



Reusable Blood Transfusion User Guide




Spike (1)

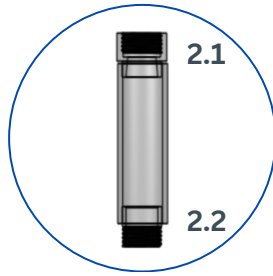


Drip Chamber
Top (2)

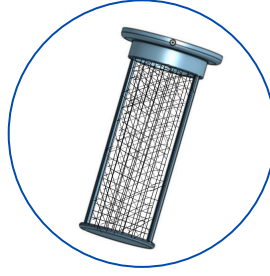

Filter Cage (3)

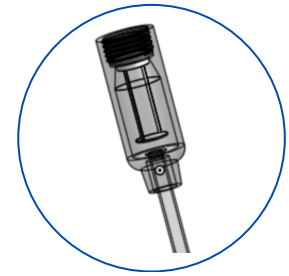


Drip Chamber
Bottom (4)

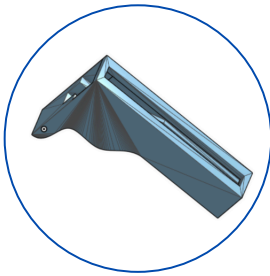

Flow Regulator (5)



Tubing (6)

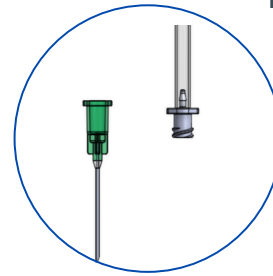

Needle Connection
(7)

Assembly

1. Obtain components listed above
2. Thread the bottom of the tubing (6) through the flow regulator (large part of flow regulator upright) (5)
3. Place filter (3) into drip chamber bottom (4)
4. Screw in the drip chamber top (2.2) to the drip chamber bottom (4)
5. Screw in the spike (1) to the drip chamber (2.1)
6. Attach the needle connection (7) to the bottom of the tubing (6)

Disassembly/ Sterilization

1. Deconstruct the transfusion set, throwing away the needle
2. Place the __ separate components into the wash prior to autoclave
3. Transfer the components into the autoclave
4. Start autoclave cycle
5. Remove components when the cycle is finished